

How to adjust contrast if you can't see anything

If contrast has been over-adjusted then here are the steps to adjust it back into the visible range. (Incidentally, in many cases if you tilt the MegaBee screen, you may still be able to see the writing on the menus, even though it is faint).

1) What MegaBee have you got?

Version 4.0 or higher



Version 3.9 or lower



The essential difference is the position of the POWER-ON button (not the colour)



Power on, pressing the respective power on button and wait for 5 seconds.

Then press the FUNCTION BUTTON once:



This will get you into the menu system.

Press:

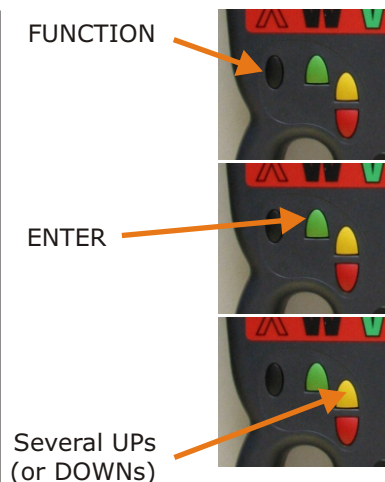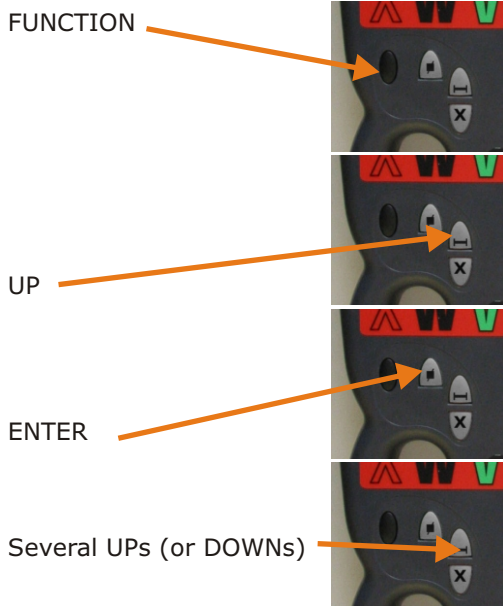

When Contrast is correct press ENTER followed by FUNCTION

